




POLIMASTER[®]
Radiation Detection Technologies

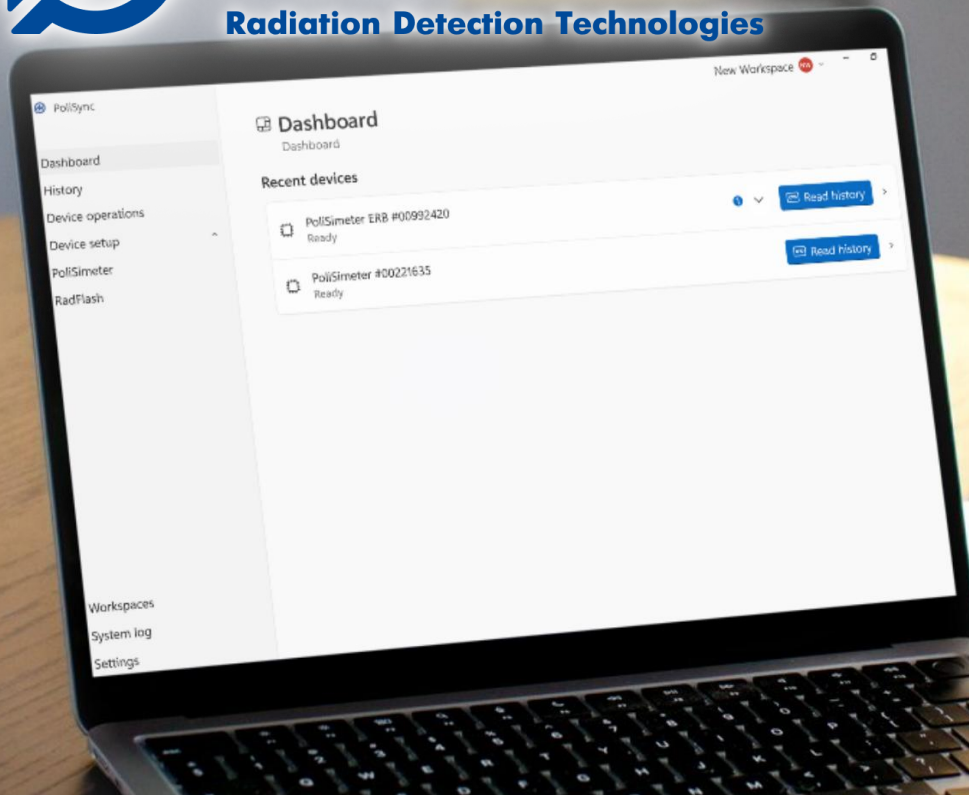

CONFIGURATION SOFTWARE

PoliSync[™]


PoliSync[™] is designed to simplify everyday work with Polimaster products by combining device setup, history readout, data storage, and operational control in a single software environment.

PoliSync is a free Windows application available in Microsoft Store. It supports **RadFlash[®]** and **PoliSimeter[™]** dosimeters. The software communicates with connected instruments via available interfaces, supports parameter management, stores device history in a local database, and provides export options for further analysis or reporting.

When a compatible device is connected via **Bluetooth** or **USB**, PoliSync automatically detects it and provides access to key workflows, including parameter configuration, history readout, configuration set application, and operation log review.

Features


- Compatible with the **RadFlash[®]** and **PoliSimeter[™]** dosimeters
- Configuration sets for faster and consistent setup across devices
- Simplified parameter configuration, including alarm thresholds
- Locks configuration to prevent unauthorized changes
- Built-in logs for traceability and operational transparency
- Role-based access for better control over protected settings
- Local database with unlimited history storage and export
- Automatic device date and time synchronization
- Workspaces for separating data and access across teams

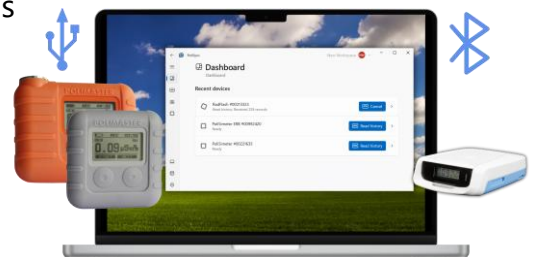

www.polimaster.com


Functionality

- Read, store, and export device history
- Adjust device parameters
- Reset accumulated dose
- Clear device history
- Create and apply configuration sets
- Create and manage multiple workspaces
- Configure access permissions for device operations

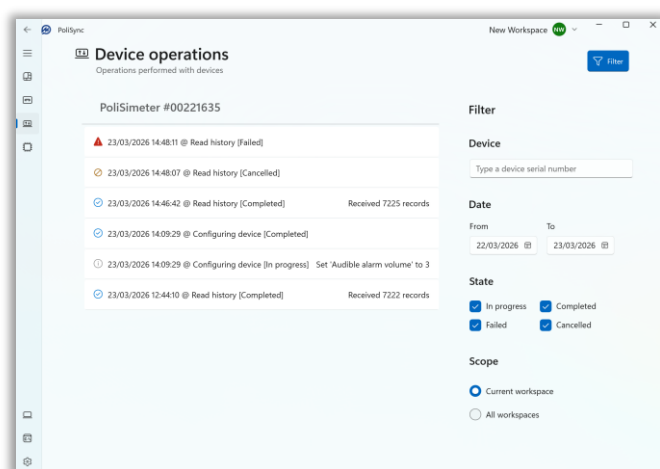


Device operations

Monitor and control all actions performed with a device in the software

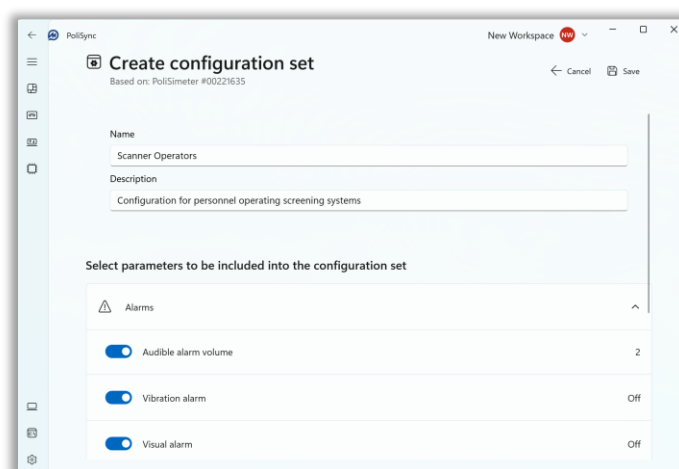

Configuration sets

Create reusable parameter templates to ensure consistent setting across multiple devices

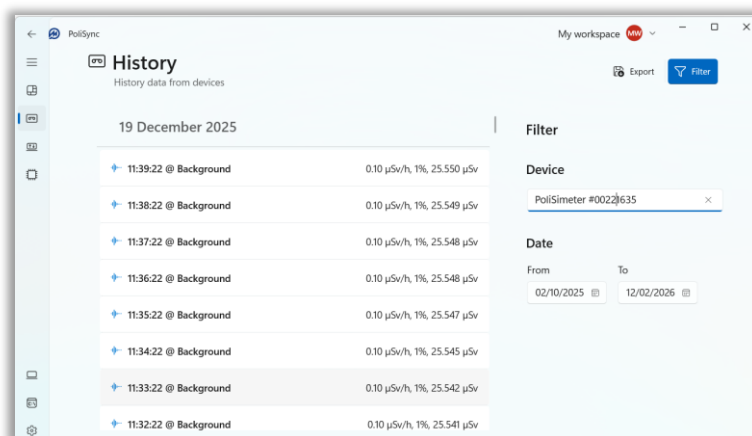

History

Read-out, store and export dosimeters history

Polimaster Inc.

45645 Willowpond Plaza,
Suite 100, Sterling, VA,
20164, USA
phone: +1 703 525 5075
fax: +1 703 525 5079
info@polimaster.us

Polimaster Europe UAB

Ezero str. 4, Didziasalio k.,
Nemėzio sen., LT-13264,
Vilnius district, Lithuania
phone: +370 5 210 2323
fax: +370 5 210 2324
info@polimaster.com

Polimaster Japan Co., Ltd.

AUBE2 5-177 Kuratsuki,
Kanazawa, Ishikawa Prefecture
920-8203 Japan
phone: + 81 076 201 8623
fax: + 81 076 201 8624
pacific@polimaster.jp